


EZ Cleaner-M

GLOSSARY OF SYMBOLS

Symbols	Meaning of Symbols
	Consult instructions for use
	Date of manufacture
	Store at 15° C - 30° C.
	Expiry date
	Manufacturer
	Keep away from sunlight
	Lot number
	Keep dry
	<i>in-vitro</i> diagnostic medical device
	Handle with care
	Direction upward

EZ Cleaner-M

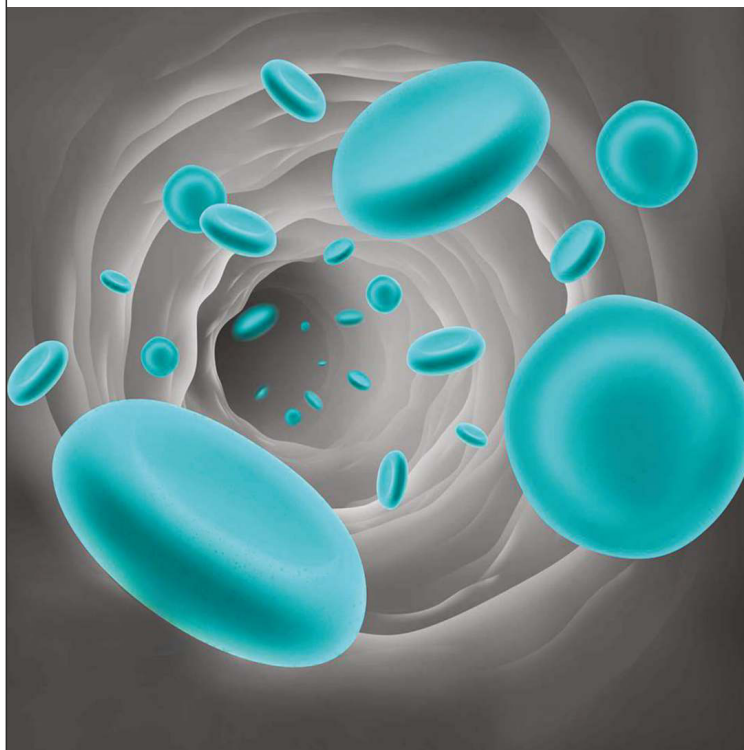

Manufactured by:

VANGUARD DIAGNOSTICS (P) LIMITED

C-123, Phase-I, Okhla Industrial Area,
New Delhi- 110 020, India.



011-41727222



info@vanguarddiagnostics.com



www.vanguarddiagnostics.com

PM300091178



EZ Cleaner-M

Pack Size : 1x100 ML



For *in-vitro* use only. Not for medicinal use.

EZ Cleaner-M

INTENDED USE

EZ Cleaner-M is used for cleaning of the tubes and baths before Shutdown in Mindray 3-Part differential Hematology Analyzers.

Please refer to the analyzer's User Manual for further information.

ACTIVE INGREDIENTS

Sodium Chloride	: 0.6%
Sodium Citrate di-hydrate	: 1.5%
Calcium Chloride di-hydrate	: 0.8%
Enzymes	: 0.02%
Dye	: 0.01%

DIRECTIONS FOR USE



- 1) Personnel installing the reagent, should be a trained laboratory professional.
- 2) Open the bottle cap of EZ Cleaner-M for carrying out the system shutdown protocol and proceed as per the instructions in the analyzer's User Manual.
- 3) Avoid contamination of the tubing and the reagent.
- 4) Do not pour and mix the remains of a reagent from a container, into another one.

EZ Cleaner-M

PRECAUTIONS

	EZ Cleaner-M, is for <i>in-vitro</i> diagnostic use only.
	EZ Cleaner-M, is environmental friendly and does not contain harmful ingredients.
	Do not inhale the reagent or pipette by mouth.
	Avoid contact with skin/eyes. In case of contact, wash off immediately with plenty of water and seek medical attention, as early as possible.
	Do not use if package is damaged.
	Keep the reagent container closed when not in use.
	Wear protective gloves when handling the reagent.
	Please refer to the MSDS associated with the reagent.
	Use Good Laboratory Practices (GLP) when handling the reagent.

EZ Cleaner-M

SPECIMEN COLLECTION AND STORAGE

Blood specimens should be collected by vein puncture in EDTA anticoagulant.

Specimens for hematological analysis may be stored at 15°- 30°C for 2 - 4 hours or at 2°- 8°C for 24 hours after collection.

REAGENT PREPARATION

This reagent is ready for use and can be connected to the analyzer directly from the container; no special reagent preparation is required.

STORAGE, STABILITY AND DISPOSAL

- Store EZ Cleaner-M between 15°and 30°C.
- The unopened reagent is stable up to the expiry date mentioned on the label when stored between 15°C - 30°C.
- Use Diluent within 06 months of first opening.
- Do not use reagent beyond the expiry date printed on the container label.
- Discard opened container after 06 months.
- Do not use reagent once frozen.
- Dispose waste product in accordance with local Government regulations.